



Value Life


INFUSION SYSTEMS




**TIVA Sets
& A Sets**
Product guide

TIVA Sets

Total intravenous anaesthesia from pump to patient

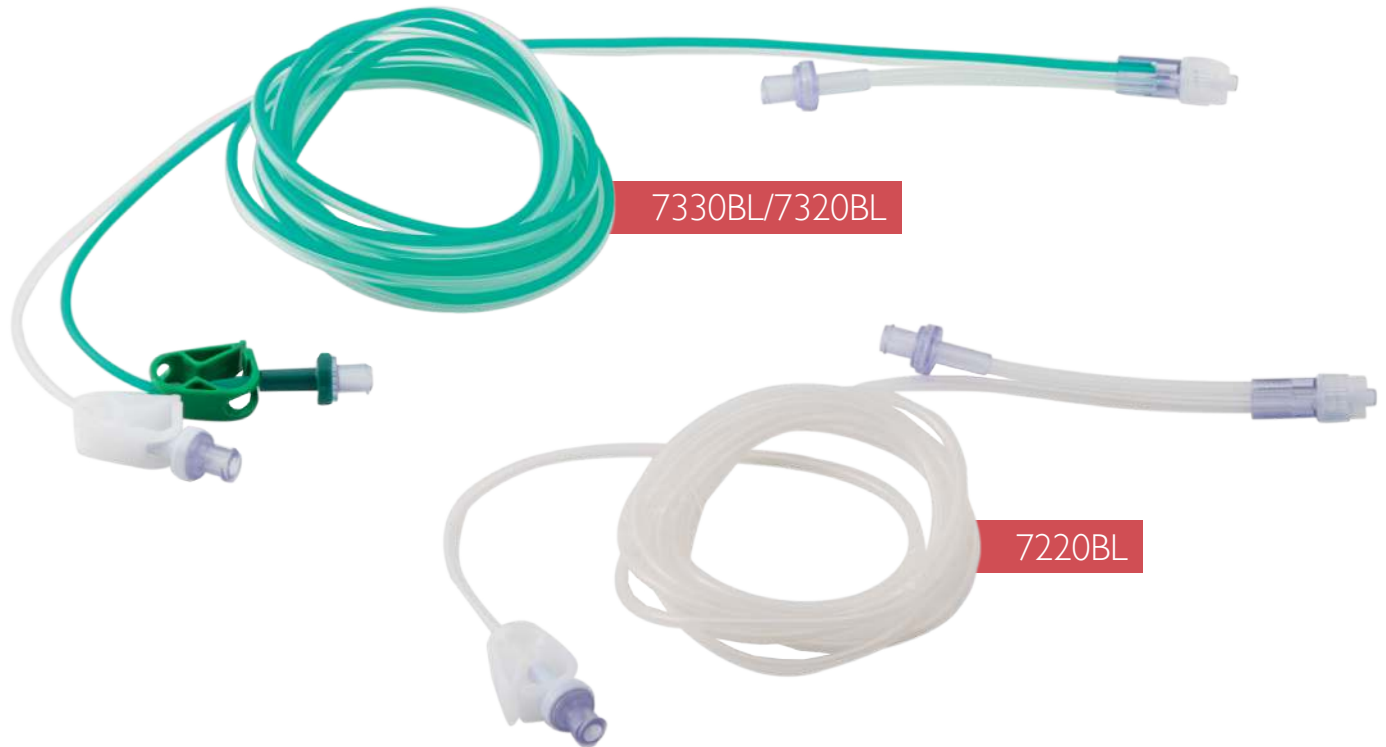


Features & benefits

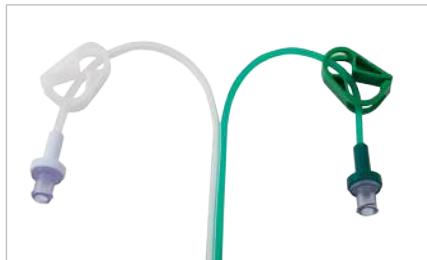

Colour-coded Siamese tubing

For clearer identification of high and low flow drug administration.

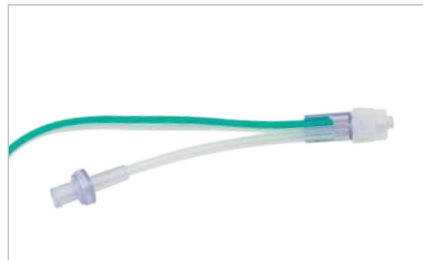

Multi-lumen from hub to tip

Keeps infusions separate, reducing the risk of drugs mixing.



Armoured tubing


Eliminates the risk of kinking, ensuring continued infusions.

Ordering codes

Vygon codes	NHSSC codes	Number of lumen	Gravity line		High flow drug line		Low flow drug line		Box qty
			Length (cm)	Priming vol. (mL)	Length (cm)	Priming vol. (mL)	Length (cm)	Priming vol. (mL)	
VAN7320BL	FSB1337	Triple	8.0	0.5	200	4.75	200	1.5	50
VAN7330BL	FSB1338	Triple	8.0	0.5	300	7.0	300	2.5	50
VAN7220BL	FSB1336	Double	8.0	0.5	200	4.75	-	-	50

A Sets

Total protection from pump to patient

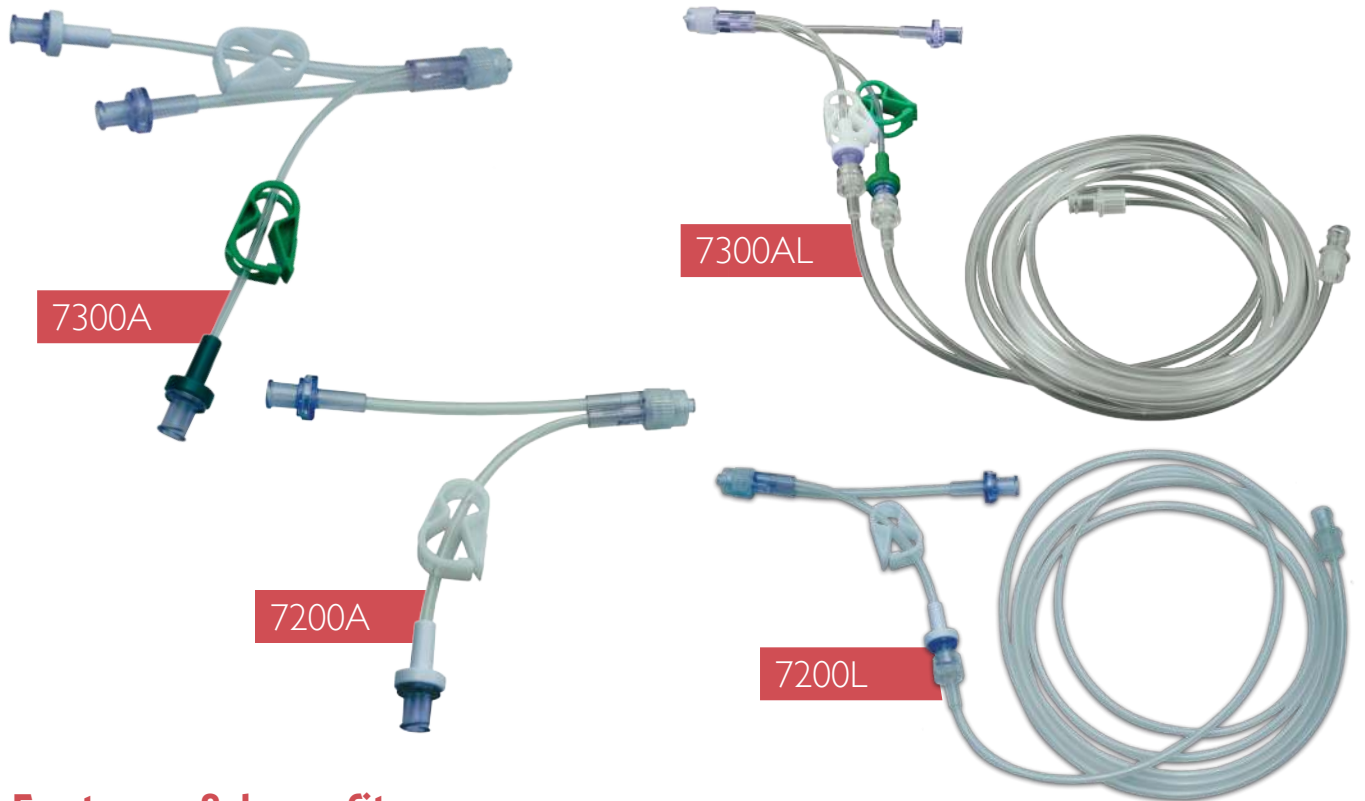

Features & benefits



Anti-syphon valve

One-way valves that require a far greater driving pressure to open, minimising the risk of downloading.



Anti-reflux valve

Minimises the risk of back-tracking and enhances patient safety.



Halkey Roberts clamp

Secure, easy-to-use clamps, which enable one-handed operation during medication changes.

Ordering codes

Vygon codes	NHSSC codes	Number of lumen	Detachable extension line	Gravity line		High flow drug line		Low flow drug line		Box qty
				Length (cm)	Priming vol. (ml)	Length (cm)	Priming vol. (mL)	Length (cm)	Priming vol. (mL)	
VAN7300A	FSB1081	Triple	No	8.0	0.5	10	0.6	12	0.3	50
VAN7300AL	FSB1082	Triple	Yes	8.0	0.5	185	5.0	187	4.7	50
VAN7200A	FSB1079	Double	No	8.0	0.5	10	0.6	-	-	50
VAN7200L	FSB1080	Double	Yes	8.0	0.5	210	5.1	-	-	50

For further information, please contact: vygon@vygon.co.uk

The specifications shown in this leaflet are for information only and are not, under any circumstances, of a contractual nature.

Vygon (UK) Ltd, The Pierre Simonet Building, V Park, Gateway North,
Latham Road, Swindon, Wiltshire SN25 4DL

Tel: 01793 748800 Fax: 01793 748899 Email: vygon@vygon.co.uk

vygon.co.uk [@vygonuk](https://twitter.com/vygonuk) [in](https://www.linkedin.com/company/vygonuk) vygonuk [yt](https://www.youtube.com/channel/UCvYgONk) vygonuk [ig](https://www.instagram.com/vygonuk) vygonuk

Copyright Vygon (UK) Ltd 2019 Content correct as of: 12/2019 Code:V00566 v1



Value Life